

Thermocycler for thermal cycles
Spectrophotometers
Colorimeters
Microtiter strip reader and water quality meter
Microplate reader and washer
pH-meters
Conductivity meter

page 238
pages 239 to 242
pages 243
page 244
page 245
pages 246 to 247
page 248



Thermocycler for thermal cycles "B960"

WITH FOUR ADAPTER MODULES FOR 0.2ML VIALS, 0.5ML VIALS AND WELLS. ADJUSTABLE TEMPERATURES FROM 0°C TO 99 °C.



Advanced technology. Peltier effect.

APPLICATIONS

It is used in molecular biology to amplify DNA via a method based in the polymerase chain reaction process.

BASIC PRINCIPLE

The equipment performs a certain quantity of thermal cycles depending on the method used and repeat them many times for a while, and at the end, the initial DNA fragment Fiber have been replicated thousands times.

For a better process output, changes between temperature levels must be made with minimum time. With thermocycler K96, the cycle temperature can be reached in seconds, even if beginning from remote positions of the last set point. These changes take place by keeping a perfect uniformity between different block points.

The system can also be programmed to produce a linear gradient of temperature widthways the block. This achieves that the highest level of productivity points of the process is optimized and located.

FEATURES


The thermocycler is composed of an inner lid system with heater and with adjustable height in order to perfectly adapt to the samples size. This prevents condensations in the upper side of the samples.

The equipment is based on a heat bomb controlled by continuous electric current and composed by several thermo electrical modules of Peltier effect, a low thermal resistance radiator and a fan convention system.


As it is integrated in the same block, this system allows increasing the process output and also fast transferring and extracting the block's temperature, going from the highest temperature level to the lowest in a minimum time.

The powerful control microprocessor allows monitoring the point where the process is at any moment, and showing it on the screen with real-time graphic images.

There is a useful elaborated software for processes organization, accessible by means of the keyboard and a high resolution LCD screen.

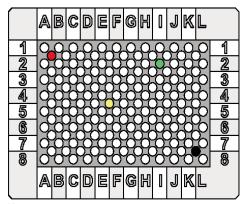

384 wells block, with red, yellow, green and black identification in diferent positions.

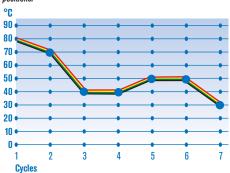

Graphic detail that shows temperature perfect uniformity in the different block wells, distinguished by different colours: red, yellow, green, black. It also shows how uniformity is maintained during all the thermal cycle periods, previously programmed in the display.



INTERCHANGEABLE BLOCKS

Each block incorporates a connector that identifies it and allows the thermocycler to recognize it. There's also an extractor handle which makes easy the block support use.



Thermal cycler door open with a heater system on the back side of the door to avoid condensations.

TECHNICAL SPECIFICATIONS

Temperature range: 0°C to 99°C.
Durability: 99min. 99sec.
Warming-up speed: 4°C/sec.
Cooling down speed: 4°C/sec.
Uniformity: At 95°C ±0.4°C.

From 20°C to 75°C ±0.2°C.

Precision: ±0.2°C.

Programmed gradient: From 2°C to 30°C according to the program.

Temperature range gradient: 30-99 °C.

Heater lid: From 70 to 115°C according

to the program.

Max. cycles number: 299.

Stored programs: Up to 1000.

Graphic display of 14.5cm, 320x240 pixels.

USB 2.0, LAN and RS232 output.

CONTROL PANEL

START switch.
Interactive graphic display.
Numeric and functions keyboard.
Height levelling control of the inner lid.

MODEL

Part No.	Height / Width / Depth (exterior) cm	Power supply	Power W	Weight Kg
5109000	25 27 38	220V/50-60Hz	780	7,8

ACCESSORIES

Interchangeable modules:

A. 96 vials of 0,2 ml. Code: **5109001**C. 96 vials of 0,2 ml + 77 vials of 0,5 ml. Code: **5109003**D. 384 wells Code: **5109004**



Visible range spectrophotometers "V-1100" and "VR-2000"

"V-1100" MODEL WITH MANUAL WAVELENGTH SETTINGS AND AUTOMATIC BLANK.

"VR-2000" MODEL WITH AUTOMATIC WAVELENGTH SETTINGS AND BLANK.



They are widely used in colleges and enterprises for general quantitatives analysis and experiments based in absorbance measurements.

COMMON FEATURES

High quality silicon photometric diode detector and 1200 lines/mm diffraction grating ensure the high quality accuracy and precision.

Digital display for an easy readout.

Automatic zero and blank (easy to use). Easy switching of transmittance, absorbance and concentration modes, just by pressing one key.

Large sample compartment, which can accommodate 5 to 100mm path length cuvettes with optional holders.

Its pre-aligned design makes it possible to change the halogen lamps by the user himself.

Optional:

Optional software based on Windows® which can expand the applications to standardcCurve and kinetics.

V-1100 MODEL

Easy of use.

Ergonomic and solid design for a continuous suitable use, for students and workers.

VR-2000 MODEL

Large LCD screen (128x64bits).

It can display a total of 50 groups of data (3 groups per screen).

It can display standard curve and kinetics curve graphic.

The system can also save the test results.

A total of 50 data groups and 10 standard curves can be saved in the RAM memory.

At most 9 standard samples can be used to establish a standard curve. The curve and the curve equation will be displayed simultaneously on the screen. The unknown concentration solutions, can be measured by the curve. If one knows the coefficient \mathbf{k} and \mathbf{b} of the formula: $\mathbf{C} = \mathbf{k} \mathbf{A} + \mathbf{b}$, one can input the value directly.


Data is stored in the memory in case of power cut.

ACCSESORIES

(see page 241).



MODEL	V-1100	VR-2000					
Part No.	4120025 4120026						
Wavelength range	325-1	000 nm					
Spectral bandwidth	5 nm	4 nm					
Optical system	Single beam, diffraction grating 1200 lines/mm						
Wavelength adjustment	Manual Automatic						
Wavelength accuracy	±2 nm	±1 nm					
Wavelength repeatability	1 nm	0,5 nm					
Photometric accuracy	±0,3% T ±0,5% T						
Photometric repeatability	±0,3% T						
Photometric range	-0,3-3 A, 0-200% T. 0-9999 Concentration						
Stray light	0,5% T	0,3% T					
Stability ±	0,004 A/h @ 500 nm	± 0,002 A/h @ 500 nm					
LCD Display	3,5 Digits	128x64 pixels					
Detector	Silico	on					
Sample compartment 10) mm path length standa	rd cuvette (100mm optional)					
Light source	Tungste	n lamp					
Output	USB, Parallel	port (printer)					
Mains supply	220 V / 50Hz AC or 110 V / 60 Hz AC						
External dimensions (HxWxD)	160 x 480 x 360 mm	180 x 470 x 370 mm					
Weight (Kg)	8	12					
SPARE Halogen lamp (visible)	V-1100 4312004	VR-2000 4312007					

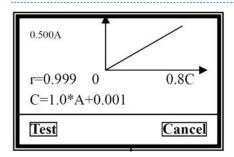

Graphic display visualized detail of a regression line for concentration calculation. (VR-2000)



Ultraviolet and visible range spectrophotometers "UV-2005" and "UV-3100"

AUTOMATIC WALVELENGTH POSITIONING AND BLANK SETTING



APPLICATIONS

Suitable for pharmaceutical and biochemical laboratories, kinetics, quantitative analysis, wavelength scanning, multiple components and DNA/protein analysis.

COMMON FEATURES

Spectrophotometers UV-2005 and UV-3100 have been developed for accurate tests; Its stray light is only of 0.05% T. They are flexible, easy-to-use and maximize value. Value is provided from every day perfomance. The independent local software provides the following functions: Basic Mode, Quantitative Test, Kinetics and System Utilities.

Data is stored in the memory in case of power cut.

The Special application software provides a complete control of the spectrophotometer from a PC, through the built-in USB port.

Automatic wavelength settings.

Halogen and deuterium lamps can be switch on/off individually to extend lamp lifetime. Its pre-aligned design makes it possible to change lamps for the user himself.

Large sample compartment for 5-100mm path length cuvettes, with optional holders. A wide range of optional accessories can also be selected.

UV-2005 MODEL

Large LCD screen (128x64bits).

It can display a total of 200 groups of data (5 groups per screen).

It can display standard curve and kinetics curve graphic.

The system can also save the test results.

A total of 200 data groups and 200 standard curves can be saved in the RAM memory.

UV-3100 MODEL

Absorbance, transmittance and concentration measurements.

It establishes or uses stored calibration equations to measure solutions of unknown concentration.

Spectrum scan of sample at any selected wavelength range with choice of scanning speed and wavelength interval.

Measurement of absorbance changing versus time with reaction rate calculation function.

Measurement at multiple wavelengths to analyse and determine the composition of the mixtures.

Calculation of concentration and DNA purity.


Note: ratio at other wavelengths can be measured.

ACCSESORIES

(see page 241).



MODELS	UV-2005	UV-3100			
Part No.	4120020	4120021			
Wavelength range	190-	1100 nm			
Spectral bandwidth	2	! nm			
Optical system	Single beam, diffraction	n grating 1200 lines/nm			
Wavelength accuracy	±0	,5 nm			
Wavelength repeatability	0,	3 nm			
Photometric accuracy	±0,	3% T			
Photometric repeatability	±0,	2% T			
Photometric range	-0,3-3 A, 0-200% T. 0-9999 Concentration				
Stray light	0,05% T @2	220 nm, 340 nm			
Stability	± 0,002 A	/h @ 500 nm			
LCD Display	Graphic(128x64)	Graphic(320x240)			
Scanning velocity	- High,Med,	Low.Max 3000 nm/min.			
Uniformity	-	± 0,002 A (200-1000nm)			
Sample compartment	10mm path length standard c	uvette (100mm optional)			
Light source	Halogen and deuteriu	m lamps (pre-aligned)			
Output	USB, parallel	USB, parallel port (printer)			
Power requirements	220 V / 50Hz AC o	r 110 v / 60 Hz A C			
External dimensions(HxWxD	180 x 470 x 370 mm	160 x 480 x 360 mm			
Weight (Kg)	14	16			

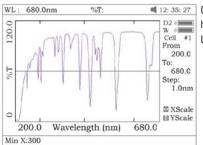

Graphic display visualized detail of holmium crystal scanning, from a UV3100 spectrophotometer.

 SPARE PARTS
 UV-2005
 UV-3100

 Halogen lamp (visible)
 4312007

 Deuterium lamp (UV)
 4312006

ACCESSORIES

Thermostatted flow Peltier cuvette holder.

Temperature range: from 15 °C. to 40 °C. in steps of 0.1 °C (ambient temperature \pm 22°C.)

Precision: ± 0,2 °C.

Sample suction volume: 1, 2, 3, 4, 5, 7.5, 10, 12.5, y 15 ml/min.

LCD Display: 128x64 pixels.

Model

Part No.	Height / Width / Depth	Height / Width / Depth	Power	Weight
	(controller) cm	(cuvette holder) cm	W	Kg
4120018	11,5 27 17	10 12 7	30	3,5

Cuvette holder

For long cuvettes, pathlength between 1 and 50 nm. Part No. 4120030 For path leength100 mm. cuvettes. Part No. 4120031

Test tubes holder

For tubes from 10 to 20 mm. Ø Part No. 4120032

Ink printer (not thermal). Paper 2 1/4" (56 mm) wide roll. Includes power transformer and interface cables. (5V-3A)

Model

Part No.	Height A	Height / Width / Depth cm			
4120117	4	16	10	1	













FEATURES

Range of polystyrene, glass and quartz cuvettes. Standard size pathlength 10mm x 45mm high. Special cuvette pathlength 40mm x 45 mm high (fig 6). Special cuvette pathlength 4mm x 45 mm high (fig 1).



MODELS

MODELO					
Standard cuvette	s				
Part No.	Figure	Material	Optical path length	Description	Presentation
5100020	9	Quartz	10 mm	Standard square, "macro"	Pack of 2 units
5100021	9	Glass	10 mm	Standard square, "macro"	Pack of 2 units
5100022	7	Polystyrene	10 mm	Standard square, "macro", disposable	Box of 100 units
5100025	8	Glass	10 mm	Flow through cell, square "micro"	Pack of 2 units

Special cuvettes					
Part No.	Figure	Material	Optical path length	Description	Presentation
5100014	6	Glass	40 mm	Rectangular	Pack of 1 unit
5100015	5	Glass	10 mm	Flow through cell. Window 5 x 10 mm	Pack of 1 unit
5100016	4	Glass	10 mm	Standard square, "micro"	Pack of 1 unit
5100017	3	Glass	10 mm	Standard square, "macro" with round lid	Pack of 1 unit
5100018	2	Quartz	10 mm	Semi-micro	Pack of 2 units
5100019	2	Glass	10 mm	Semi-micro	Pack of 2 units
5100023	7	Polystyrene	10 mm	Square Semi-Micro, disposable	Box of 100 units
5100024	1	Glass	5 mm	Square Semi-Micro	Pack of 1 unit







NECTA: Double Beam Spectrophotometer "UV-2300 II"

BAND WIDTH 1.5 nm. TOTAL STABILITY DOUBLE BEAM OPTICS. AUTOMATIC WAVELENGTH CORRECTION CALIBRATION. USB PORT FOR DATA STORAGE. WIDE RANGE OF ACCESSORIES.





APPLICATIONS

Research, chemistry, biotechnology, general spectroscopy analysis applications, environmental applications.

FEATURES

Monochromator high resolution optics that eliminate any optical aberrations, monochromator "Seya-Namioka", manufacturer of exclusive beam technologies and diffraction gratings in Japan.

Bandwidth 1.5 nm in accordance with European pharmacopoeia recommendations. (the relation between the maximum and minimum absorbence in Toluene and Hexane at 0.02% (V/V) should be more than 1.5T).

Several modes of operation including spectral scanning, time base scanning, multi wavelength determinations, peak and trough detection, etc.

Fast spectra scan displayed on a screen covering the whole spectral range: quick scan 3600 nm/minute, range 190 to 1100 nm

Validation function for GLP/GMP:

This function maintains and assures the optimum working parameters of the instrument. Parameter such as the wavelength precision and noise are monitored.

Memory facilities of analytical results:

The analytical parameters and results can be stored in the "flash" memory, connected by the USB. The stored information can be stored as text and can be transferred to a computer for reporting using MS WORD/EXCELL.

A DNA/RNA function is fitted to quantify the ratio at 260/280 nm.

Controllable from a computer with application specific optional software: "UV-Analyst Spectrum", (see accessories).



SPECIFICATIONS

Optical system: Double beam optics. Wavelength range 190 nm to 1100 nm.

Band pass: 1.5 nm

ABS range: -2.000 to 3.000 A or 0 to 300% T.

Stray light: less than 0.05% (220 nm Nal, 340 nm NaNO₂). Scan speed: 10, 100, 200, 400, 800, 1200, 2400 and 3600 nm/minute.

Wavelength accuracy: ±0.3 nm.

Photometric accuracy: ± 0.002 A from 0 to 0.5 A.

± 0.004 A from 0.5 to 1.0 A.

Baseline stability: 0.003A/hr (500 nm after 2 hours of use.).

Noise level: 0.003 A (at 500 nm).

Light source: Deuterium D2 and Halogen lamps.

Built in screen: LCD back light of 190 x 138 mm.

Connections: RS232 and parallel port.





Constant temperature 6 place cell changer, 10mm cuvettes, ideal for Kinetics.

Part No. 5110029



USB port. "Flash" memory not included.

SPARES

Tungsten halogen lamp. Part No. **5110021**

Deuterium lamp (UV).Part No. 5110022

AU	U	Eg	90	K	13
_					

Part No. **5110023** Flow cell 10 mm path length.

Part No. 5110024 Constant temperature single cell holder for 10mm flow cell.

Part No. 5110025 Micro cell holder for 10mm (50 µl) cell.
Part No. 5110026 Long path length cell holder for 100mm cells.

Part No. 5110027 5 Position cell changer for 10mm cuvettes.
Part No. 5110028 Constant temperature cell holder for 10mm

Part No. 5110028 Constant temperature cell holder for 10mm cuvettes
Part No. 5110029 6 place constant temperature cell changer for 10mm cuvettes.
Part No. 5110033 Software UV-Analyst Spectrum. Simple and logical, enhanced

software for multiple applications and results and manipulation of

data, calculation such as DNA/RNA ratios.



Semi Automatic clinical analyser "Photometer S-2000"

LOW PROFILE FUNCTIONAL CLINICAL ANALYSER.

APPLICATIONS

Clinical laboratories, Clinical Biochemistry, Haematology, E.I.A. Electrolytes and lones.

FEATURES

47 most commonly used Pre-programmed clinical methods: cholesterol, glucose, blood iron, LHD, AST, UREA etc.

Storage for up to 134 methods and up to 1000 stored results.

6 analytical fitting routines: final point, bio-chromatic, regression fit, multi-standard curve, best fit and kinetics.

Graphic display LCD (640 x 240 pixels, 256 colours).

Lamp life save function.

Graphic display.

Quality control of analysis function.

Built in printer.

Control by "roller ball mouse" or "external mouse."

Temperature conditions: 15 °C to 32 °C: Maximum humidity 85%

RS232 out put interface.

SPECIFICATIONS

Solid state photodiode detector.

Bichromatic optical system with 7 filters, 340, 405, 500, 546, 578, 620 and 670 nm.

Band pass: 10 nm. Flow cell: 30 ul.

Measuring range: from 0.000 to 2.500 ABS.

Resolution: display 0.001 ABS, calculations to 0.0001ABS.

Thermostatted Peltier cuvette holder: 25, 30 and 37 °C Precision ±0.1 °C.

Tungsten halogen lamp of 6V 10 W.

Comes complete with protective cover, spare lamp and two rolls of printer paper.



MODEL

Part No.	Height / Width / Depth cm		/ Depth	Power rating	Power W	Weight Kg	
4120010	14	45	33	110-220V/50-60Hz	150	10	

ACCESSORY

Metallic thermostat dry block "BIO-BATH"

Small size with a fixed temperature of 37 $^{\rm o}{\rm C}$ with check thermometer, 8 tube positions for 75x13 VAC tubes.

Note: Accelerates the extraction of plasma for faster results and waiting time for the patient.

MODEL

Part No.	Height / Width / Depht			Stability	Power	Weight
	(exterior) cm			°C	W	Kg
7001561	8.8	9.6	7	0.5	8	0.5

SPARE. Halogen lamp 6 V 10 W. Part No. 4312015



Digital colorimeter "Clormic"

MICROPROCESSOR CONTROLLED.
AUTOMATIC ZERO ABSORBANCE AND 100% TRANSMITTANCE.
ALPHANUMERIC 20 CHARACTER 2 LINE L.C.D. DISPLAY.

APPLICATIONS

Chemistry Laboratory, quality control, environmental control.

FEATURES

Wavelength range: 400 to 800 nm, by using special filters with a 30 nm band pass.

12 position Filter disc wheel selection.

Standard filters: 420, 440, 490, 520, 550, 580, 620 and 680 nm.

Expanded Absorption range: -0.3 to 3.5 O.D. in real time.

Transmission: 0 to 100 T %.
Photometric accuracy: >1 %.
Photometric precision: ±1 %.
Photometric stability: 0.004 A/hr.
Light source: Long life tungsten lamp.

Detector: Solid state.

Sample chamber: 10 mm cuvettes or test tubes of 12 and 16 mm \emptyset .

Minimum volume: 1 ml.

Display: Alphanumeric LCD display of 2 lines of 20 characters.

Calculation functions: Transmission T %. Absorbence, Concentration by factor or standard concentrations.

Calibration: Self adjusting by software. **RS-232 interface.**

CONTROL PANEL

ON/OFF switch.

Interactive LCD display. Numeric and function keypad.



MODEL

Part No.	Built	Height / Width / Depth	Power	Weight
	in printer	cm	W	Kg
4120009	NO	11 18 28	10	4.5

SPARE

Lamp of 6 V / 6 mm. Part No. **4512009**



Microtiter Strip Reader "Reader M-2000"

STORAGE MEMORY FOR UP TO 59 METHODES INCLUDING: FACTORS, STANDARDS OR CALIBRATION CURVES WITH UP TO 10 STANDARDS. ADDITIONAL CALCULATION FEATURES.

APPLICATIONS

Clinical diagnostics, E.I.A., Veterinary, quality control.

FEATURES

Spectral range: 400 to 700 nm, using interference filters of 10 nm band pass.

Supplied complete with a filter of 450 nm (other wave-

lengths available consult for options.)
Measuring range: -0.2 to 3.5 0.D.
Photometric accuracy: >1%.
Photometric precision: ±1%.
Light source: Long life LED.

Detector: solid state. 59 method memories. Curve point editing feature.

Software: ABS. Concentration and polynomial curve fitting with up to 10 standards.

RS-232 output to connect a printer or to a PC.
Interface RS-232.



MODEL

Part No.	Height / Width / Depth		Power	Weight	
	cm		W	Kg	
4120013	11	18	28	7	2

ACCESSORY

Ink printer (not thermal). Paper 2 1/4" (56 mm) wide roll. Interface RS232. Includes power transformer and interface cables. (5V-3A)

Model

Part No.	Height / Width / Depth cm	Weight Kg
4120113	4 16 10	1



SPARE

Lamp of 4 mm Ø 5 V. Part No. **4313030**



Multi-parameter water quality meter "MP-2006"

LARGE TOUCH SCREEN. TO MEASURE PH, ORP, CONDUCTIVITY, TDS, SALINITY, DISSOLVED OXYGEN, IONS & TEMPERATURE

FEATURES

5.7 inches touch screen.

Capacity to store up to 2200 sets of data. GLP method. Possibility of storing, printing and deleting measuring data. USB and RS-232 ports for PC connection.

Measuring modes: continuous, timing and balance measuring modes.



	Range	Accuracy	Resolution
pH:	-2 to +19,999	±0,002 FS	0,001
ORP:	±1999,99 mV	±0,03 FS	0,01 mV
Conductivity:	0,001 µS to +199,9mS	±0,5 FS	0,001µS
TDS:	0,000 mg/L to 19.99g/L	±0,5 FS	0,001 mg/L
Salinity:	(0,0 to 8,00)%	±0,1%	0,1%
Dissolved oxygen:	0,00 to 19.99 mg/L	±0,5 mg/L	0,01 mg/L
lons:	0 to 19990(ug/L, mol/L)	±0,5%	±1
Ta:	-0,5 to +135	±0,5 FS	0,1



MODEL					
Code	Height	/ Width cm	/ Depth	Power	Weight Kg
4100700	20	20	7	OOON/EO COU-	4

ACCESSORIES

Standard electrodes of pH, conductivity and temperature. See page 247.

Special electrodes:

oheeiai eid	outouto.								
Code	Electrode type	Working temperature	Body material	Measure mm	Code	Electrode type	Working temperature	Body material	Measure mm
4120183	Platinum	0-50⁰	Glass	Ø 12 x 120	4120189	Calcium	5-60⁰	Glass	Ø 12 x 120
4120184	Platinum-2	9 0-50°	Glass	Ø 12 x 120	4120190	Nitrate	5-60⁰	Glass	Ø 12 x 120
4120185	Antimony	0-50⁰	ABS	Ø 12 x 120	4120191	Copper	5-60⁰	Glass	Ø 12 x 120
4120186	Tungsten	0-50⁰	Glass	Ø 12 x 120	4120192	Chloride	5-60⁰	Glass	Ø 12 x 120
4120187	Silver	0-50⁰	ABS	Ø 12 x 120	4120193	Oxygen	5-40⁰	ABS	Ø 12 x 120
					4120194	Reference	5-55⁰	Glass	Ø 12 x 120



Microplate reader 2100-C

TOTALLY AUTOMATIC TOUCH SCREEN CONTROL.



APPLICATIONS

Clinical diagnostics. Quality control in food analysis, Haematology.

FEATURES

Bi-chromatic optical system with 4 wavelength selection:

 $405,\,450,\,492,\,630$ nm, (other filters from 405 to 700 nm available). Memory for up to 500 programs, 10000 sample results and 1000 patients

Patient data software, laboratory personnel and analysis reports.

100 typical analysis methods pre programmed.

Capable of performing up to 12 different tests per plate.

 ${\it Calculations: ABS, Cut-Off, \ Curve \ fitting: regression, \ linear, \ logarithment \ logar$

mic, exponential and potential.

Measuring absorbance range: 0-2.5A

Displayed absorbance range: 0-3.5A

Accuracy:±1.0% or ±0.007 A.

Precision: ±0.5% or ±0.005 A.

Linearity: r > 0.995.

Resolution:0.001 ABS (display); 0.0001 ABS (calculated).

Measuring time: Continuous mode < 5 s, step by step mode < 15 s. Output RS-232 for download to a computer or a printer.

Parallel printer interface.

Programmable plate vibrating for mixing with control of time and speed.

CONTROL PANEL

LCD display (320 x 240 pixels) of 5.7 inches.

Operating system WINDOWS CE.

Operation by touch screen or optional "Mouse" control (mouse not in-

cluded).



MODEL	Part No.		Width / Depth erior) cm	Voltage	Weight Kg
2100-C	5109999	19	45 33	AC 110-250 V / 50-60 Hz	10

ACCESSORIES

Software for data transfer to computer. Part No. 5109997



COMECTA Microplate washer 2600-C

COMPLETELY AUTOMATIC. EASY TO USE. 50 WASH PROGRAMS.



APPLICATIONS

Clinical laboratories.

Blood banks.

Quality control food laboratories.

FEATURES

Wash head of 8 x 1. (12 x 1 optional).

Plate wash by 12 x 8 strips or 8 wells/strip.

Plate clean for flat, round and V bottom plates.

Automatic wash and rinse processes.

Reservoirs: 2 litre, wash solution.

2 litre, waste, with liquid level sensor.

Memory for up to 50 wash programs.

Software controll of head's depth adjust and horitzontal position.

Wash cycle times from 1 sec. to 2 hrs. Wash dosing range: from 50 to 2000ul.

Dosing resolution: 5% CV (based on 350ul of distilled water)

Residual volume (after aspiration):

< 1 µl (Cells, bottom U and V well plates.).

< 5 μ l (Cells, bottom of flat well plate.)

CONTROL PANEL

LCD graphic display of 90 x 50mm Membrane key pad of 9 keys. Emergency stop button.



MODEL	Part No.	Height / Width / Depth (exterior) cm	Voltage	Weight Kg
2600-C	5110000	14 45 33	AC 110-250 V / 50-60 Hz	12



Digital pH-meter "pH-2003" PORTABLE METER.

FEATURES

Manual temperature compensation from 0 to 60 °C with sample temperature key pad entry (does not read temperature).

Automatic buffer recognition of two of those three buffers: pH4, pH7 and pH9.

Electrode connector: BNC type.

Includes: carry case and pH electrode.

SPECIFICATIONS

	Range	Precision	Resolution
рН:	0 to 14	±0.03	0.01
mV:	±1400	±2	1

Impedance: 5 x 10¹¹ Ohm.

MODEL

Part No.	Height / Width / Depth cm	Power source	Weight Kg
4120300	17 7.5 3	2 x 1.5V batteries	0.5



Digital pH meter "pH-2004"

PORTABLE METER WITH IP65 PROTECTION AGAINST WATER AND DUST.

FEATURES

Automatic temperature compensation and read out from 0 to 60 °C with temperature sensor Rt-2252. Results memory for up to 250 readings.

Protection against water and dust IP65 (Impermeable to water splashes from any direction.) Automatic buffer recognition and calibration for 1 or 2 buffers, recognition of up to 5 selectable buffers.

RS-232 output to connect a printer or a computer. Electrode connector: BNC type. Includes: carry case, pH electrode and temperature sensor Rt-2252.

Direct connection to printer Part No. 4120113 (see page 247).

Universal arm support not included. Part No. 1001552 (see page 248).

SPECIFICATIONS

	Range	Precision	Resolution
pH:	0 a 14	±0.02	0.01
mV:	±1800	±2	1
Ta:	0 to 60	±0.5	0.1
Impedance:	1 x 10 ¹² Ohm.		

MODEL

Part No.	. 3	idth / Depth m	Power source	Weight Kg
4120400	20 1	0 4.5	4 x 1.5V AA Batteries	0.5



Digital pH meter "pH-2005"

BENCH TOP MODEL WITH LARGE SCREEN.

Automatic or manual temperature compensation from 0 to 99.9 °C Automatic buffer recognition and calibration of 2 buffers, recognition of up to 3 buffers: pH 4, pH 7 and pH9.21.

RS-232 output to connect a printer or a computer.

Electrode connection BNC.

Temperature sensor connection Rt-2252.

Includes electrode's support stand.

SPECIFICATIONS

	Range	Precision	Resolution
pH:	0 to 14	±0.02	0.01
mV:	±1999	±1	1
Ta (°C):	0 to 99.9	±0.5	0.1
Imnedance:	1 x 10 ¹² Ohm.		



See "Accessories for pH meters" for selection of electrodes (Page 242).



ACCESSORY

pH probe and additional temperature sensor, gel filled. Temperature sensor Rt 2252. With IP65 protective connectors.

Part No. 4120401



ACCESSORY

Portable magnetic stirrer, battery operated and submersible. Part No. 7001731 for "pH-2003" and "pH-2004" (see page 27).







MODEL

Part No.	Height / Width / Depth cm	Weight Kg
4120500	30 20 7.2	1.5

See accessories for electrode, temperature sensor and buffer solutions. Not compatible with printer Part No. 4120113 (see page 247).



Digital pH meter "pH-2006"

BENCH TOP, LARGE GRAPHIC TOUCH SCREEN CONTROL.

FEATURES

Ergonomic touch screen control functions. Includes support arm.

Automatic temperature compensation from -5 to 105 °C.

Up to 5 buffer calibration with 10 buffer recognition.

Stores up to 200 results.

RS-232 output to connect a printer or a computer.

Electrode connector: BNC type. Rt-2252 Temperature port. Selectable resolution range.

Direct interface for printer Part No. 4120113 (see accessories).

SPECIFICATIONS

	Range	Precision	Resolution
pH:	-2 to 18	±0.01	0.001
mV:	±1999.9	±0.1	0.01
T ^a (°C):	-5 to 105.0	±0.3	0.1

Impedance: 1 x 10¹² Ohm.

See accessories for electrodes, temperature probes, buffers and printer.



	_	_	÷	÷
IVI	Ш	Ш	H	ı
	_	_	-	•

Part No.	Height / Width / Depth cm	Weight Kg	
4120600	29 20 7	1	

ACCESSORIES for pH meters.





Universal arm support.

With holder for up to 4 sensors or electrodes. Very stable heavy base.

Ideal for working at various heights.

Part No. 1001552

General laboratory use.



μ	П	Electiones
F	ig.	Part No.

Fig.	Part No.	Range pH	Sensor T ^a	Temperature Range	Electrolyte	Dimensions mm	Material Body
1	4120102	0 - 14	No	0 - 80°	Liquid	Ø 12 x 120	Glass
2	4120125	0 - 14	Yes	0 - 60°	Gel	Ø 12 x 120	Ероху
3	4120101	0 - 14	No	0 - 60°	Gel	Ø 12 x 120	Ероху
4	4120114	0 - 10	No	0 - 60°	Gel	Ø 12 x 120	Ероху
5	4120104	0 - 14	No	0 - 60°	Gel	Ø 12 x 120	Glass
6	4120103	0 - 14	No	0 - 80°	Gel	Ø 8 x 60	Glass
7	4120126	0 - 14	Yes	0 - 80°	Liquid	Ø 12 x 120	Glass

у	Protective
у	Protective

Application

cover. General laboratory and external use. cover. General laboratory and outside use.

Surface analysis, Skin, meat, moist semisolids, culture media. Pointed sensor. Suitable for food, semi-solids, fruit meat etc.

For measuring in tubes. Liquid solutions in general.

Fast response, liquids in general.

Temperature sensors

Fig.	Part No.	Range	Туре	Dimensions mm	Material
8	4120121	0-100	Rt-2252	Ø 5 x 120	Stainless steel
9	4120120	0-100	Rt-2252	Ø 12 x 120	Glass

Ion selective electrodes

IOII GOIGGE	ion octoonive electroace								
Part No.	lon	Range	Limits	Range °C	Range pH				
4120168	Chloride Cl ⁻	1-3 x 10 ⁻⁶	35.500 - 1	5-50	1-12				
4120174	Nitrate NO⁻₃	1-7 x 10 ⁻⁶	62.000 - 0.4	0-50	2-11				
4120177	Potassium K⁺	1-10-6	39.000 - 0.04	0-50	1-9				

pH Buffer solutions.

4120107 Flask 250 ml Buffer ±0.02 pH 4.00 at 20 °C; 4.01 at 25 °C Flask 250 ml Buffer ±0.02 pH 7.02 at 20 °C; 7.001 at 25 °C 4120108 Flask 250 ml Buffer ±0.02 pH 9.23 at 20 °C; 9.18 at 25 °C 4120127 Frasco 250 ml KCI 3M solution (Electrolyte for pH electrodes) 4120109

MAGNETIC STIRRER "AGIMICRO" Features

Maximum stir volume: 2 litres. Adjustable stir speed from 40 to 1400 r.p.m.

Excellent chemical resistance.

Location hole at the back for a retort rod. (Accessory)

Ink printer (Not thermal).

interface cables. (5V-3A)

For use with pH and conductivity meters . Paper 2 1/4" (56 mm) wide roll. Interface RS232. Includes power transformer and



Control panel

Mains switch with luminous on indicator. Analogue control of speed in r.p.m.

Model Part No. Medidas Power source

Kg 230 V 1 W 7001638 Ø 12 x 5 Supplied complete with a 8 Ø x 32 mm P.T.F.E. coated stir bar.



ACCESSORY Support bar and electrode clamp. Part No. 7001639

ŀ	Model			
	Part No.	Height / Width / Depth cm	Power source	Weight Kg
l	4120113	4 16 10	DC 5 V, 3 A	1

Weight



Conductivity meter "CD-2004"

PORTABLE EQUIPMENT WITH 1P65 PROTECTION AGAINST WATER AND DUST.

FEATURES

Measures conductivity, TDS and salinity.

Automatic range switching.

Automatic temperature compensation.

Calibration using standards for conductivity, TDS and salinity.

RS-232 output to connect a printer or a computer.

Direct printer connection, see part no. 4120113.

Temperature coefficient programable.

Conductivity cell constant adjustment.

Memory for 250 results.

Ambient protection IP65 (protection against water splashes from any direction).

Includes Carry case, conductivity cell K= 1 and temperature sensor Rt-2252.

Conductivity cell with epoxy body. Temperature sensor Rt-2252.

Conductivity cell and temperature sensor.

IP65 protection connectors. Part No. 4120411

SPECIFICATIONS

Reading ranges: (automatic) Conductivity

0.00 to 19.99 μ S/cm 0.00 to 10.00 mg/L

20.0 to 199.9 µS/cm 10.00 to 100.0 mg/L 200 to 1999 µS/cm 100.0 to 1000 mg/L 2.00 to 19.99 mS/cm 1.000 to 10.00 g/L 20 to 199.9 mS/cm 10.00 to 19.99 g/L (Cell constant K=10) (Cell constant K=5 or 10)

Salinity range: 0.00 to 8.00 Temperature range: 0.0 a 40.0 °C

Precision:

Conductivity: 1.5% base scale

0.20% Salinity: Temperature: 0.3 °C



MODEL

Part No.	Width / Height / Depth			h Power	Weight
	cm			source	Kg
412041	0 21	10	4.5	4 1.5 V batteries	0.5



ACCESSORY

Conductivity meter "CD-2005"

BENCH TOP METER WITH LARGE SCREEN.



Measures conductivity and TDS.

Automatic or manual setting temperature compensation with fixed coefficient of 2%.

Calibration with standards of conductivity and TDS. RS-232 output to connect a printer or a computer.

Cell constant adjustment.

Temperature sensor connection Rt-2252. Temperature measuring range: 0 to 60 °C.

SPECIFICATIONS

Conductivity ranges:

(Manual range change)

0 to 20.00 µS/cm

0 to 200.0 μ S/cm 0 to 2000 µS/cm

0 to 10.00 mS/cm

Precision: 1.5% base scale + 1 digit. TDS ranges: 0 to 1000 mg/L.



30 20 7.2

Accessories to complement the meter, conductivity cells, temperature sensors and standards (see accessories). Not compatible with the printer 4120113.

4120510

Includes sensor stand.



Conductivity cell

Glass body. Working temperature range 0-100 °C

l	Fig.	Part No.	K	Sensor I" ("C)	Dimensions mm	Applications
	1	4120220	1	Yes	Ø 12 x 120	General laboratory use up to 2000 µs
	3	41202220.	1	Yes	Ø 12 x 120	Pure water up to 20 µS
	4	4120223 1	0	Yes	Ø 12 x 120	High conductivity up to 200 ms

Temperature sensors

Fig.	Part No.	Range	Туре	Dimensions mm	Material
5	4120121	0-100	Rt-2252	Ø 5 x 120	Stainless steel
6	4120129	0-100	Rt-2252	Ø 12 x 120	Glass

Conductivity standards.

4120160 Flask 250 ml Standard ±0.05 of 1278 mS/cm at 20 °C; 1413 mS/cm at 25 °C 4120161 Flask 250 ml Standard ±0.025 of 4.915 mS/cm at 20 °C; 5.446 mS/cm at 25 °C 4120162 Flask 250 ml Standard ± 0.06 of 11.67 mS/cm at 20 °C; 12.88 mS/cm at 25 °C

MAGNETIC STIRRER "AGIMICRO" Features

Maximum stir volume: 2 litres.

Adjustable stir speed from 40 to 1400 r.p.m.

Excellent chemical resistance.

Location hole at the back for a retort rod. (Accessory)

Control panel

Mains switch with luminous on indicator. Analogue control of speed in r.p.m.

Model

Part No.	Medidas cm	Power source	Weight Kg
7001638	Ø 12 x 5	230 V 1 W	0.5

Supplied complete with a 8 Ø x 32 mm P.T.F.E. coated stir bar.



0.5

ACCESSORY Support bar and electrode clamp. Part No. 7001639



Universal arm support.

With holder for up to 4 sensors or electrodes. Very stable heavy base. Ideal for working at varying heights.

Part No. 1001552